

Please contact Omni Customer Service with your sizing details so we can make sure your shipment is sent correctly.

1-802-497-2253 - OPT 1

or

customer care@urincare.com

SIZING CHART

Small	25mm
Medium	29mm
Intermediate	32mm
Large	36mm
X-Large	41mm

STYLE CHART




TYPE 1 (FULL) - Provides full hydrocolloid adhesive coverage along the entire sheath for a secure, reliable fit. Designed to help maintain placement during daily activities while supporting skin integrity. Recommended for extended wear of up to 24-48 hours, as clinically appropriate.



TYPE 2 (SHORT)
Features full hydrocolloid adhesive coverage on a shorter sheath design, making it ideal for users with a shorter penile length or for shorter wear durations. Commonly recommended for nighttime use or for patients who meet the minimum anatomical length requirements while still requiring a secure, comfortable fit



TYPE 3 (MID BAND)
Features a full-length sheath with a reduced hydrocolloid adhesive band, providing secure placement while minimizing skin contact with adhesive. Designed for non-active or non-ambulatory patients who require reliable urine management with less adhesive coverage. Recommended for users with limited mobility where maximum adhesive retention is not necessary.

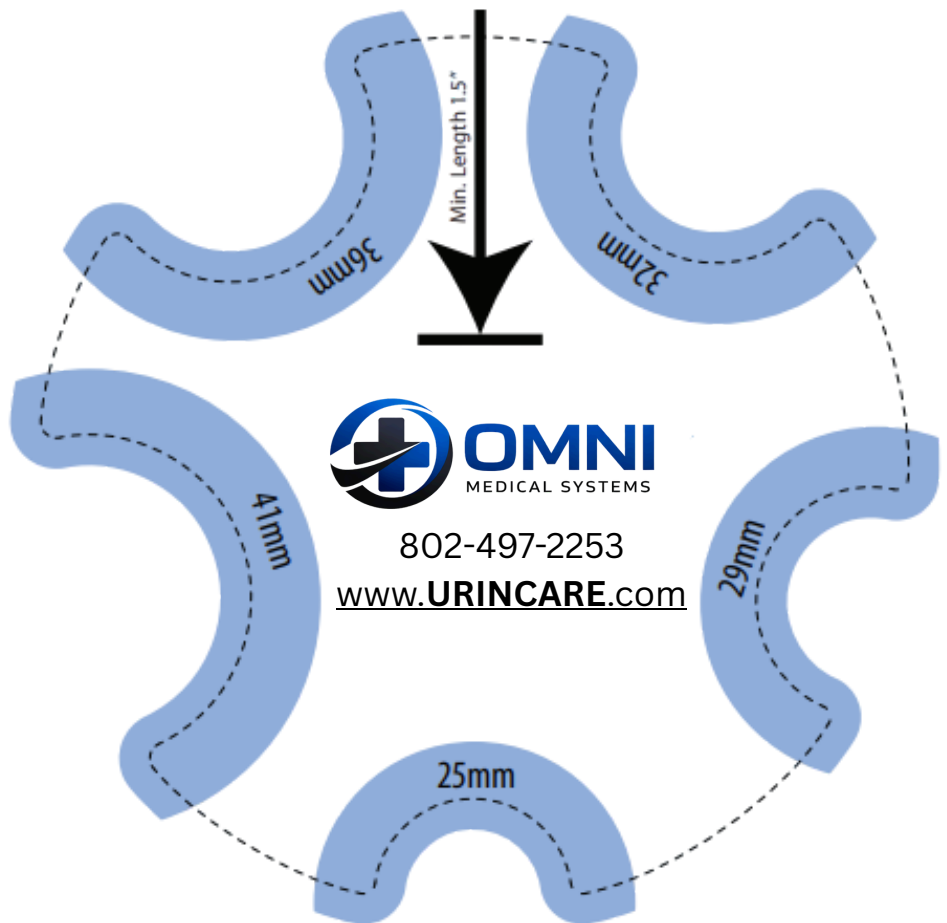

802-497-2253

www.URINCARE.com



Cut along the dotted line to remove the URINCARE FLO™ Sizer.
Use the sizer to determine the appropriate size for proper fit and performance.

URINCARE FLO™ SIZING INSTRUCTIONS

Important: This sizer is intended for reference purposes only. Please contact Omni Medical Systems to obtain an official URINCARE FLO™ Sizer for accurate sizing and product selection.

1. Carefully place the sizer around the base of the penis and adjust it to find the most comfortable and secure fit. This measurement will determine the appropriate sheath diameter.
2. Caution: Always remove the penis from the sizer before rotating or adjusting the sizing tool.
3. Once the diameter has been determined, position the sizer at the base of the penis and verify that the penile length extends beyond the minimum required length of 1.5 inches.
4. Record the diameter measurement and consult the sizing chart to select the appropriate URINCARE FLO™ size.

Note: Proper sizing is essential for optimal fit, comfort, adhesion, and performance of the URINCARE FLO™ system.

Bad Fit



Ideal Fit

