

APPLYING THE URINCARE FLO™ CONDOM

(STEP BY STEP GUIDE)

Prepare your hands

Wash, rinse, and dry your hands completely before handling the condom.

Wash and clean the area

Wash the penis gently with soap and water. Rinse thoroughly and pat dry.

Check the skin

Inspect the penis for any redness, broken skin, or irritation. Do not apply if the skin is not healthy—seek medical advice if unsure.

Trim excess hair (if needed)

Carefully remove hair at the base of the penis to help ensure a secure fit and prevent pulling.

Position correctly

Hold the penis gently at a 90-degree angle away from the body.

Place the device

Position the URINCARE FLO™ at the tip of the penis, leaving about a half-inch space between the tip of the penis and the FLO connector.

Roll into place

Slowly and gently roll the condom portion back over the penis, ensuring it fits smoothly and securely.




Connect the tubing

Attach the hose to the FLO connector. You will hear a soft “click” or snap when it is securely locked in place.

To disconnect

When removing, gently squeeze the hose connector while holding the FLO connector steady, then pull them apart carefully.

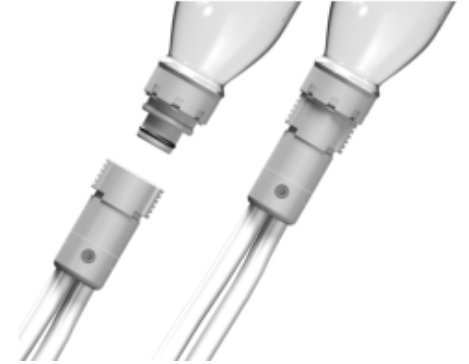

Connect the Hose to the Collection Bag

Attach the free end of the hose to the port on control unit on the collection bag. Ensure the connection is secure before use.

Automatic Urine Collection

The URINCARE FLO™ system operates automatically. When urine is detected by the integrated sensors, the control unit activates the pump, gently drawing urine away from the skin and into the collection bag. This helps keep the user dry and promotes skin integrity.

Intermittent Catheter Access

The hose can be easily disconnected from the FLO™ Connector when intermittent catheterization is required. This allows catheter access without removing the FLO™ condom, helping maintain proper placement and reducing the need for reapplication.

