

WHY CHOOSE URINCARE™?

- ✓ No Adhesives Required
 - Ideal for men with sensitive skin
 - Eliminates adhesive-related irritation and skin breakdown
- ✓ Actively Keeps You Dry
 - Sensor-driven technology automatically detects urine
 - Quiet pump transfers urine away from the body into a collection bag
- ✓ Promotes Skin Health
 - Reduces prolonged skin exposure to moisture
 - Helps minimize skin irritation and complications associated with incontinence
- ✓ Discreet & Comfortable
 - Worn comfortably under clothing
 - Portable, lightweight, and battery-powered
- ✓ Freedom to Stay Active
 - Can be worn while sitting, standing, walking, or semi-reclined
 - Supports an active and independent lifestyle

WHO MAY BENEFIT?

The URINCARE™ Male **CUP** System may be appropriate for:

- ✓ Veterans living with urinary incontinence
- ✓ Post-prostatectomy patients
- ✓ Individuals with Multiple Sclerosis (MS), ALS, Parkinson's Disease, SCI, or TBI
- ✓ Men with sensitivity to adhesives
- ✓ Patients who have experienced skin breakdown from traditional products

WHAT'S INCLUDED IN THE CUP STARTER KIT?

- 1 - Sensor-Driven Control Unit
- 2 - Rechargeable Batteries
- 1 - Charger with power cord adaptor
- 3 - Male CUP™ Body Interface
- 3 - Adjustable Comfort Harness
- 3 - Collection Bags

VETERAN SUPPORT

**Covered by the Veterans Administration
for Eligible Veterans**

Talk with your VA Provider or Prosthetics Department to learn if URINCARE™ is right for you.



U.S. Department
of Veterans Affairs

CONTACT US

OMNI MEDICAL SYSTEMS

808 HERCULES DRIVE
COLCHESTER, VT 05446, USA



PHONE: 1-802-497-2253

WWW.URINCARE.COM



URINCARE
— NONINVASIVE BLADDER MANAGEMENT SYSTEM —
★ FOR VETERANS ★




**MALE CUP™
QUICK START GUIDE**





CONTROL UNIT

Automatically activates when sensors in the CUP detect urine, gently pumping urine away from the skin to help keep the user dry.

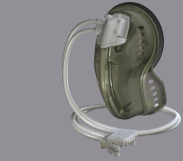

MALE CUP

The CUP is ideal for men with sensitive skin, skin irritation, or existing wound care concerns. Its non-adhesive design minimizes skin contact while helping keep users dry, comfortable, and protected. Designed for everyday use, the CUP can be worn while sitting, standing, walking, or semi-reclined.



COLLECTION BAG

Dual-chamber collection system helps reduce skin contact and irritation. The leak-proof Hydro-Block™ air filter helps keep users dry and control odor. Available in waist, leg, and bedside options for added convenience.



MALE HARNESS

The URINCARE™ Harness securely holds the CUP in place without the use of adhesives, helping ensure a comfortable, reliable fit. Soft, adjustable Velcro® straps allow users to customize the fit for optimal comfort and proper CUP positioning throughout daily activities.



STEP 1: CHARGING THE BATTERY

1. Insert the rechargeable **battery** into the **charger**.
2. Connect the **charger** to the wall adapter and plug it into an outlet.
3. Check the status light:
 - ● Solid Green = Fully Charged
 - ● Flashing Yellow = Charging
 - ● Flashing Red = Service Required

Note: Fully charge the battery before first use. Initial charging may take up to 2.5 hours.

STEP 2: CONNECT THE BATTERY

1. Connect the charged **battery** to the **Control Unit** (slide within the grooves), until it locks into place.

STEP 3: SELECT THE COLLECTION BAG

1. Choose which **collection bag** that best fits your needs
 - 600mL or 1000mL **Leg Bag** (worn on the thigh or calf)
 - 1000mL **Waist Bag**
 - 2000mL **Bedside Bag**



STEP 4: MALE CUP & UNDERGARMENT PLACEMENT

1. Insert the **Male Cup** into the **undergarment** with the foam ring facing the body.
2. Place the penis into the opening and ensure it is fully positioned inside the **cup**.
3. Secure the **Undergarment** by fastening the adjustable Velcro® straps on both sides. Adjust as needed to achieve a secure, comfortable fit and ensure proper positioning of the Male Cup.

STEP 5: COLLECTION BAG PLACEMENT FOR USE

1. Secure the **Collection Bag** using the option that best fits your needs:
 - Leg Bag: Attach to the leg using the provided Velcro® straps.
 - Waist Bag: Secure around the waist using the provided buckle and adjust the straps for comfort and a secure fit.
 - Bedside Bag: Position the bedside bag in a convenient location as needed.

STEP 6: CONNECTING THE SYSTEM FOR USE

1. Slide the **Control Unit** (pump) into the grooves on the bag clip - on the **Collection Bag**, until it clicks and locks into place. Once secured, it cannot be removed without pressing the black release button.
2. Attach the **Male Cup Connector** to the **Control Device** (already connected to the **Collection Bag**) by sliding it into the guide grooves until fully seated. The **Control Device** will beep twice and the battery status lights will flash to confirm proper connection.



THAT'S IT! YOU'RE IN CONTROL WITH URINCARE™!

ENJOY GREATER INDEPENDENCE AND CONFIDENCE WITHOUT THE CONSTANT WORRY OF URINARY INCONTINENCE.

The flashing lights on the Control Device indicate if the unit is fully charged, partially charged or drained.

LED Status:

- **Green:** Full Battery or 100-81%
- **Green/Yellow:** 80-61%
- **Yellow:** Medium Battery 60-41%
- **Yellow/Red:** 40-21%
- **Red:** Battery Low or 20-0%

The Male CUP™ can be worn while sitting, standing, walking, or in a semi-reclined position.

