

WHY CHOOSE URINCARE™?

- ✓ Designed for Active Lifestyles – Wear while sitting, standing, walking, reclining, or sleeping
- ✓ 24/7 Urine Management – Continuous, around-the-clock protection
- ✓ Hydrocolloid Adhesive Comfort – Secure, skin-friendly fit for extended wear
- ✓ Patented FLO™ Technology – Promotes airflow while actively removing urine
- ✓ Keeps Skin Dry – Helps reduce moisture buildup and skin irritation
- ✓ Smart Sensor Activation – Automatically detects urine and activates the system
- ✓ Intermittent Catheter Friendly – Remains in place during catheterization
- ✓ Kink-Resistant Tubing – Maintains consistent flow and helps prevent backflow
- ✓ Reliable Leak Protection – Supports skin integrity and dependable performance
- ✓ Greater Freedom & Independence – Flexible, discreet, and designed to keep you active and confident

WHO MAY BENEFIT?

The URINCARE™ Male FLO System may be appropriate for:

- ✓ Veterans with urinary incontinence
- ✓ Post-prostatectomy and prostate cancer survivors
- ✓ Overactive bladder (OAB) and urge incontinence
- ✓ Spinal Cord Injury (SCI), MS, Parkinson's, ALS, or TBI
- ✓ Neurogenic bladder patients
- ✓ Veterans who leak between intermittent catheterizations
- ✓ Active or ambulatory veterans seeking an alternative to briefs and pads
- ✓ Veterans requiring day and nighttime urine management
- ✓ Side, stomach, or back sleepers
- ✓ Veterans seeking a discreet, hands-free bladder management solution
- ✓ Veterans looking to reduce caregiver burden and improve independence

WHAT'S INCLUDED IN THE CUP STARTER KIT?

- 1 - Sensor-Driven Control Unit
- 2 - Rechargeable Batteries
- 1 - Charger with power cord adaptor
- 15 - URINCARE FLO™ Condom Catheters (up to 48 hours wear time)
- 3 - Hoses with Sensors
- 3 - Collection Bags

VETERAN SUPPORT

Covered by the Veterans Administration for Eligible Veterans

Talk with your VA Provider or Prosthetics Department to learn if URINCARE™ is right for you.

VA



U.S. Department
of Veterans Affairs

CONTACT US

OMNI MEDICAL SYSTEMS
808 HERCULES DRIVE
COLCHESTER, VT 05446, USA



PHONE: 1-802-497-2253

WWW.URINCARE.COM



CONTROL UNIT

Automatically activates when sensors in the CUP detect urine, gently pumping urine away from the skin to help keep the user dry.

FLO (CONDOM CATH)

The URINCARE FLO™ features a patented Hydro-Block™ odor-resistant air filter that promotes continuous airflow, helping prevent vacuum pressure while supporting smooth, unrestricted urine drainage. Applied with medical-grade hydrocolloid adhesive provides a secure, skin-friendly fit for extended wear.

COLLECTION BAG

Dual-chamber collection system helps reduce skin contact and irritation. The leak-proof Hydro-Block™ air filter helps keep users dry and control odor. Available in waist, leg, and bedside options for added convenience.

MALE HARNESS

Patented, kink-resistant hose with integrated sensors that detect urine and automatically activate the system to move urine from the FLO™ condom into the collection bag. Helps maintain consistent flow, keep the user dry, support skin health, and reduce the risk of backflow.



STEP 1: CHARGING THE BATTERY

1. Insert the rechargeable **battery** into the **charger**.
2. Connect the **charger** to the wall adapter and plug it into an outlet.
3. Check the status light:
 - Solid Green = Fully Charged
 - Flashing Yellow = Charging
 - Flashing Red = Service Required

Note: Fully charge the battery before first use. Initial charging may take up to 2.5 hours.

STEP 2: CONNECT THE BATTERY

1. Connect the charged **battery** to the **Control Unit** (slide within the grooves), until it locks into place.

STEP 3: SELECT THE COLLECTION BAG

1. Choose which **collection bag** that best fits your needs
 - 600mL or 1000mL **Leg Bag** (worn on the thigh or calf)
 - 1000mL **Waist Bag**
 - 2000mL **Bedside Bag**



STEP 4: APPLY THE URINCARE™ FLO (Condom)

1. Wash Hands - Wash, rinse, and dry your hands thoroughly.
2. Clean the Area - Wash the penis with soap and water. Rinse and dry completely.
3. Inspect Skin - Check for redness, irritation, or broken skin. Do not apply if skin is damaged.
4. Trim Hair (If Needed) - Trim hair at the base of the penis to improve adhesion and comfort.
5. Position the Penis - Hold the penis at a 90° angle away from the body.
6. Apply the **URINCARE FLO™** - Place the **FLO™** at the tip of the penis, leaving approximately ½ inch of space between the tip and the FLO™ connector.
7. Roll into Place - Gently roll the **FLO™** down the penis, ensuring a smooth, secure fit with no wrinkles.




STEP 5: COLLECTION BAG PLACEMENT FOR USE

1. Secure the Collection Bag using the option that best fits your needs:
 - **Leg Bag:** Attach to the leg using the provided Velcro® straps.
 - **Waist Bag:** Secure around the waist using the provided buckle and adjust the straps for comfort and a secure fit.
 - **Bedside Bag:** Position the bedside bag in a convenient location as needed.

STEP 6: CONNECTING THE TUBING TO THE FLO™

1. Attach the **hose** to the **FLO connector** (while the condom has already been applied to the penis). You will hear a soft "click" or snap when it is securely locked in place.

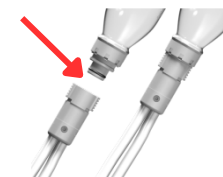

STEP 7: CONNECTING THE SYSTEM FOR USE

1. Slide the **Control Unit** (pump) into the grooves on the bag clip - on the **Collection Bag**, until it clicks and locks into place. Once secured, it cannot be removed without pressing the black release button.
2. Connect the **Hose** to the **Collection Bag** - Attach the free end of the **hose** to the port on **control unit** on the **collection bag**. Ensure the connection is secure before use.



The flashing lights on the Control Device indicate if the unit is fully charged, partially charged or drained.

LED Status:

Green: Full Battery or 100-81%

Green/Yellow: 80-61%

Yellow: Medium Battery 60-41%

Yellow/Red: 40-21%

Red: Battery Low or 20-0%



INTERMITTENT CATHETER ACCESS

The hose can be easily disconnected from the FLO™ Connector when intermittent catheterization is required. This allows catheter access without removing the FLO™ condom, helping maintain proper placement and reducing the need for reapplication.

THAT'S IT! YOU'RE IN CONTROL WITH URINCARE™!

ENJOY GREATER INDEPENDENCE AND CONFIDENCE WITHOUT THE CONSTANT WORRY OF URINARY INCONTINENCE.